

**FARO S.p.A.** founded in 1948 by Osvaldo Favonio, designs and manufactures equipment for dental units, dental practices and dental labs. Manufacturing, design, research & development, all take place in our over 5,000 square metre plant in Ornago (MB), Italy.

Over the years, the company has become a world leader in this specific sector: a role gained by commitment to technology, research and design and an ability to understand the real needs of end users.

Utmost quality is assured by the internal management of the whole manufacturing cycle and a qualified and efficient pre and after-sales service. FARO experience best expresses the pride of a product made entirely in Italy.

The **six product lines** that manufacture and market offer a **comprehensive choice of solutions** for equipping dental practices and laboratories.



#### **ILLUMINATION**

LED AND HALOGEN

includes several models of lamp, made with LED and HALOGEN technology;



#### **APPLICATION SETS**

AND ACCESSORIES

essential for mounting the lights and their accessories to the dental unit and within the practice;



#### **STERILIZATION**

AND PACKAGING

includes class B autoclaves, distillers and sealing machines;



#### **INSTRUMENTS**

FOR DENTAL CABINETS

turbines, syringes, contra-angles and handpieces;



#### **COMPONENTS**

FOR DENTAL UNITS

a range of accessories and parts dedicated to dental units;



#### **INSTRUMENTS**

FOR LABORATORIES

motors and handpieces;

We can also customize our products to meet specific requests. For a customization project please contact our Customer Service.



EXPERIENCE AND INNOVATION  
SINCE 1948

#### **FARO S.p.A.**

via Faro, 15 - 20876 Ornago (MB) - Italy Tel. +39 039.68781 - Fax +39 039.6010540  
[www.faro.it](http://www.faro.it) - [export@faro.it](mailto:export@faro.it)

#### **FARO FRANCE**

Za Tgv Coriolis - 71210 Montchanin - France Tel. +33 385.779680 - Fax +33 385.779688  
[www.farofrance.com](http://www.farofrance.com) - [farofrance@farofrance.com](mailto:farofrance@farofrance.com)

#### **FARO DEUTSCHLAND GMBH**

Gewerbepark Heideckhof Heideckstr. 179 D-47805 Krefeld - Germany Tel. +49 2151.936921 - Fax +49 2151.936933  
[www.faro.it](http://www.faro.it) - [info@farodeutschland.de](mailto:info@farodeutschland.de)

Certified Company



ISO 9001:2008 ISO 13485:2003

INFO LINE:



W9200017\_REV1\_ENG

LED  
TECH



**SIDÈREA** **NEW**  
**DENTAL**  
**LED ROOM LIGHT**



EXPERIENCE AND INNOVATION  
SINCE 1948

# SIDÉREA DENTAL LED ROOM LIGHT

LED  
TECH

24 MONTH WARRANTY



THE FIRST **FARO** ROOM LIGHT DESIGNED FOR DENTAL PRACTICES.

Always in search of the best lighting system for dental practices, FARO has designed **SIDÉREA**: the multi-directional room light that ensures an harmonious and homogeneous illumination of the workplace.




**Elegant, minimal design: pleasing to the eye.**

Sophisticated and functional, **SIDÉREA** creates a pleasantly elegant atmosphere to work in. The white aluminium surgical frame with chromed metal finishes matches well with the FARO family feeling and is perfectly in harmony with ALYA, the operating lamp.

**Homogeneous illumination: the perfect balance.**

**SIDÉREA** is the multi-directional room light that provides homogeneous and well distributed illumination, **leaving no area in the shadow.**

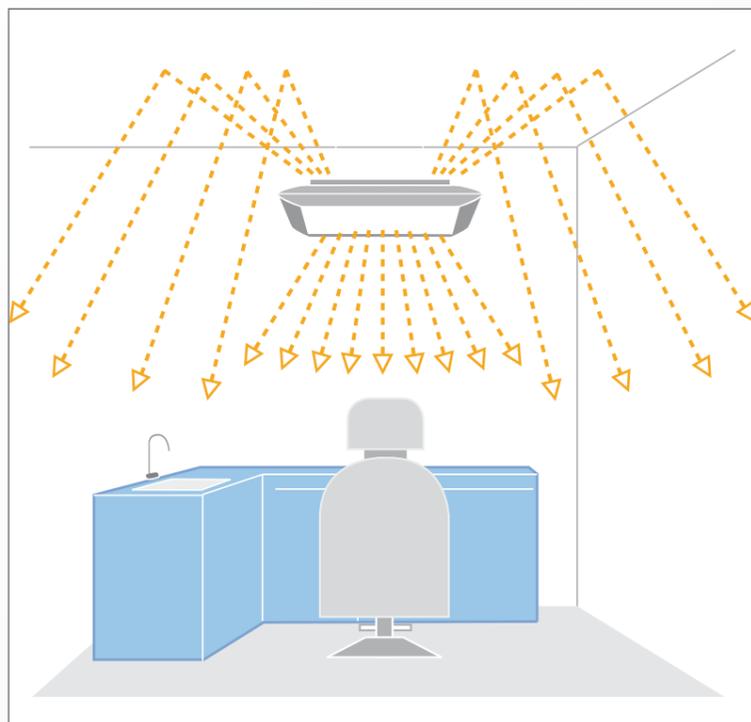

**No glare: it loves your eyes.**

**SIDÉREA** guarantees the visual comfort of patients and operators thanks to an **UGR index below 10.**



#### WHAT IS UGR?

UGR (Unified Glare Rating) is the quantitative measure of the glare of a light source within a room and ranges between 10 (no glare) and 30 (considerable glare).



**Relaxing or performing light: the colour you need.**

Thanks to the **different colour temperatures**, **SIDÉREA** allows to create the best working environment to suit all needs.



**How much light? As much as you need.**

The illumination control allows to **adapt the light to the needs of each operator**, depending on the lighting conditions of the room.



#### Remote control

**5 lighting options** pre-set by FARO, which can be reconfigured by users to meet their specific needs.

#### SIDÉREA App

Compatible with iOS and Android operating systems, to control the room light directly from any **smartphone or tablet.**

## FARO LIGHTING SYSTEM: THREE INTERACTING LIGHTS

**SIDÉREA** is designed to interact with the dental lamp **ALYA** and the pre-operative light **THEIATECH**: the three light sources can work in sync, creating the right balance of light in the room.

A complete and integrated lighting system that allows operators to move their visual focus from one working area to the other without stressing their eyes.



## TECHNICAL SPECIFICATIONS

GENERAL SPECIFICATIONS*	
Power	300 W max
Dimensions	1.595 x 635 mm
Weight	18,5 Kg
Colour Temperature	Indirect: 4.000 K / Direct: 6.500 K
Colour Rendering Index	> 95
Lighting value	up to 3.000 LUX (at 1,4 m distance)
UGR	< 10

\* (The technical data specified refer to typical values subject to tolerance)